



The Surgical Panel SP1-24 with an integrated monitor as well as build-in IT and video management components is a digital image viewing system, designed for use in the operating room. Integrated into the hospital's IT structure, it is ideal for work performed by OR nursing staff. The 24-inch monitor offers maximum user-friendliness and allows for perfect viewing of all information with ergonomic head and body position.

- Built-in 24" monitor (touchscreen option)
- Workstation for OR nursing staff
- Space-saving wall mounting
- High ease of service and maintenance
- Flexible choice of material, configuration and accessories
- Completely sealed within the OR and disinfectible



Surgical Panel SP1-24

Individually customized housing

The housing dimensions can be customized, the housing material can be rendered in steel, stainless steel, glass, or a combination of these materials. The color selection is usually from the RAL range but other options are also available.

Various monitors

Much freedom when selecting number, size and type of the integrated monitors. Partially they can be configured as touchscreens.

Flexible IT and video management

The modular interior design allows customer-specific configuration and uncomplicated servicing.

Dust and water proof

The Surgical Panel is completely sealed within the OR, thus it can easily be cleaned and disinfected. With intelligent heat distribution over the entire front, it does not even require a ventilation system.

Varied accessories

Digital watches, speakers, cameras, on-air lights, barcode readers, as well as docking stations for peripheral devices such as audio players, and data storage devices can be integrated. Furthermore, a tailor-made interface integration is possible.

Easy Maintenance and Service

The trunk principle allows the device to be opened without complicated disassembly, because of its unique locking system no screw connections have to be opened. This makes the maintenance of installed components easy and efficient.


Different Mounting Options

Due to the device's construction, surface-mounted or flush-mounted installation is possible, depending on assembly options on site.

Specifications

Model Variations		SP1-24: without touchscreen SP1-24T: with touchscreen
Housing	Material	Powder-coated stainless steel
	Front	Combination of stainless steel and anti-reflective safety glass (ESG)
	Color	Available in every RAL color
	Mounting Options	Suitable for flush mounting or surface mounting
Panel	Size	24"
	Diagonal	61 cm (16:10 aspect ratio)
	Panel Type	Color TFT LCD Panel (IPS) with DICOM Preset
	Native Resolution	1920 x 1200 (2,3 MP)
	Display Colors	16.7 million colors
	Pixel Pitch	0.270 mm x 0.270mm
	Brightness (typical)	300 cd/m ²
	Contrast Ratio (typical)	1000:1
	Viewing Angles (H / V, maximum)	178°, 178°
	Touchscreen	Projected capacitive, multi-touch (SP1-24T only)
	Response Time (typical)	5 ms (Midtone)
Backlight	LED	
Signal Inputs	DisplayPort x 1, HDMI x 1, DVI-D x 1, DSub x 1	
IT Components		INTEL® Core™ i5 processor - 3.x GHz, INTEL® Chipset, 6 MB L3 Cache, 8 GB DDR3 RAM (extendable), SATA HDD or SATA SSD, LAN RJ 45 - 10/100/1000 Base TX
Operating System		MS Windows 7 prof., Windows 10
Keyboard / Mouse		Anti-microbial medical silicone keyboard with integrated mouse storage (foldable) silicone mouse (visual 5 keys, laser mouse /IP65)
Inputs		USB 2.0 / 3.0 x 2 (front side)
Push Buttons		Main switch x 1 (front side) System switch x 1 (front side)
Maximum Power Consumption		157.7 W
Power Requirements		100 - 240 V / 50 - 60 Hz / 4 - 2 A
Thermal Load		48.75 kJ/h
Net Weight		appr. 36 kg
Certifications & Standards (Please contact EIZO for the latest information)		IEC/EN60601-1 (3rd edition), RoHS, WEEE

Dimensions (unit: mm)

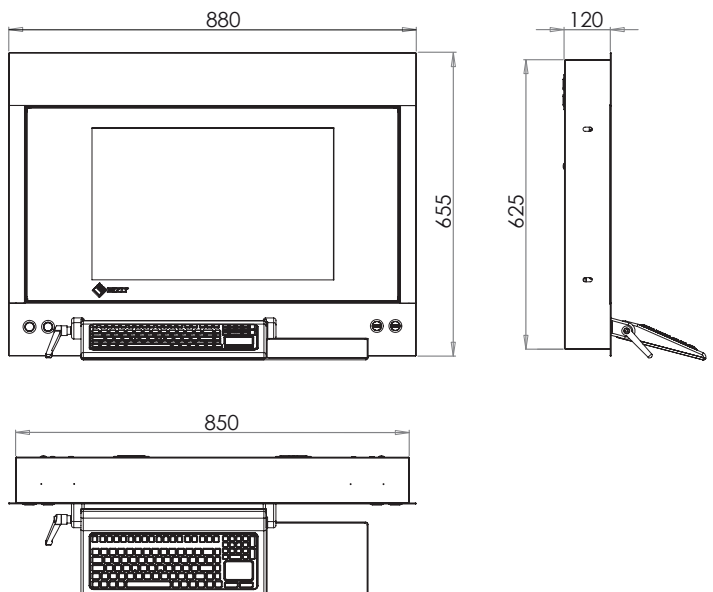

You will find your EIZO contact partner at:
www.eizo-or.com/contact

All product names are trademarks or registered trademarks of their respective companies. Curator is a trademark and EIZO is a registered trademarks of EIZO Corporation. Specifications are subject to change without notice.

Copyright © 2017 EIZO GmbH. All rights reserved. (1704)